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PL

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INSTRUCTION MANUAL EN


WARMING BELLY STRAP
COZYWARM GREY

**PLEASE READ THIS USER MANUAL BEFORE USING THE DEVICE!
FOLLOW ALL SAFETY INSTRUCTIONS TO AVOID DAMAGE
CAUSED BY IMPROPER USE! KEEP THE USER MANUAL FOR FUTURE REFERENCE.**

SAFETY INSTRUCTIONS

1. Do not use while sleeping.
2. Please avoid using, in the same position for a long time and pay attention to confirm the skin condition, if there is any abnormal immediately removed.
3. The ban on near the fire, and all kinds of corrosive liquid used in adverse environment.
4. Forbidden to use gravity disease and heart disease.
5. Skin wound is still unhealed or it is prohibited to use in patients with skin abnormalities.
6. Gout onset period or not completely disabled survivors have suffered gout. Pregnant women, the cervical tumor patients disabled.

PRODUCT DESCRIPTION

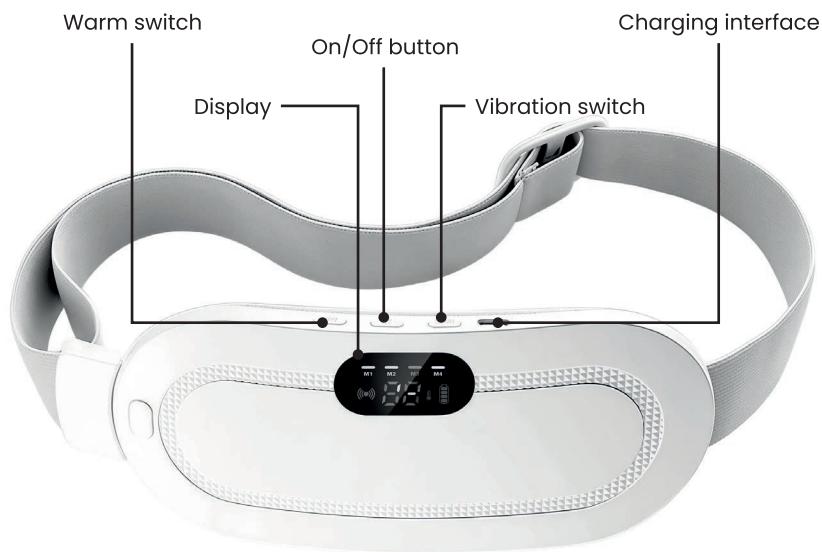USER INSTRUCTIONS

Machines using graphene heating principle, three temperature, **a 50° daily warm palace, second gear 55° period nursing three gears 60° dysmenorrhea**, can also through the vibration.

Touch the vibration switch: double-click the touch up fast, short touch switch modes, long touch to turn it off.

Warm touch switch: double-click the touch up fast, short touch switch temperature, long touch to turn it off.

MAINTENANCE AND CLEANING

If the heated menstrual belt gets dirty, gently wipe the outer fabric with a slightly damp cloth. Before cleaning, make sure the device is switched off and disconnected from the power supply. Do not immerse the belt in water or use harsh chemicals.

TECHNICAL DATA

RATED POWER	3W
VOLTAGE	DC 5V
CHARGING	5V
WARM GEAR	50/55/60°C
MATERIAL	ABS + deer velvet fabrics
SIZE	185 x 85 x 25 mm

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

Malis B. Machoński Sp. K. declares that device described in this manual meets all the requirements of the:

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