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PL

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INSTRUCTION MANUAL EN


BODY MASSAGER
SKINBOOST


PLEASE READ THIS USER MANUAL BEFORE USING THE DEVICE! FOLLOW ALL SAFETY INSTRUCTIONS TO AVOID DAMAGE CAUSED BY IMPROPER USE! KEEP THE USER MANUAL FOR FUTURE REFERENCE.

SAFETY INSTRUCTIONS

- Do not use the device if you are pregnant or suffering from heart problems or any physical injuries, and consult your physician if you are unsure.
- Only use the device one hour after meals.
- Keep away from children.
- Do not use the device on the same area of body for more than 15 minutes at a time.
- Do not use the device in bathroom or any humid places.
- If the device work disorder, please immediately cut power off, and contact the product after-service center.
- Do not use the device when you feel unwell.
- Confirm the device stop working when exchange the message heads.
- If the power wire is damaged, please ask the after-sales maintenance technician to change it.
- Do not immerse in water. Keep away from liquids or heat sources.
- Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than 30 minutes.
- Unplug the charger from device battery after charging or prior to use.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

PRODUCT DESCRIPTION

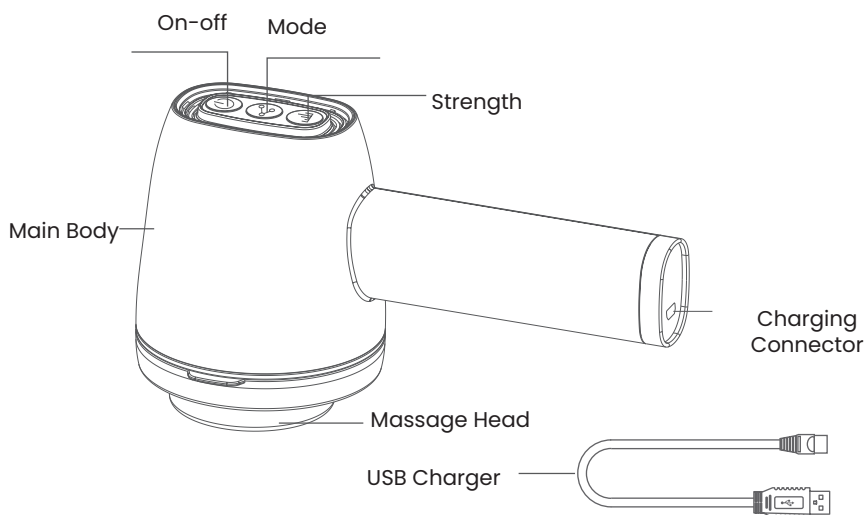
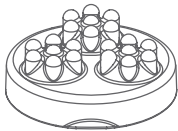
Multi-finger head


Wave head

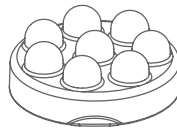Ball head

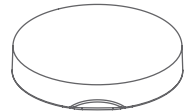Flat head

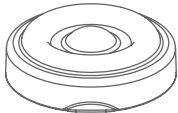Annulus head

Petal head

Brush Head

Frosted head

USER INSTRUCTIONS

1. Fully charge the product if low battery.
2. Press the Power Button, turn on the device, it is working.
3. Press the Mode Button to change the mode according your requirement, 5 modes for choice.
4. Press the Strength Button to adjust the speed according your requirement, 6 levels of speed for choice.
5. For the first time to use better with a lower massage speed, gradually increase the speed after adaption.
6. To keep the massager in best using condition, it is suggested to use more than 15 minutes each time, it will auto off when working continuously for 15 minutes.
7. After usage, turn off the power.

Charge Instructions

1. Please fully charge the product before first time using or after not use for a long time.
2. Power off the product before charge it by provided USB cable.
3. Power light is flashing during charging and transfer to long lighting after fully charged.
4. Unplug the USB cable from charging connector after fully charged.

MAINTENANCE / CLEANING AND STORAGE

1. Use soft cloth or gentle cleanser to clean machine. Massage heads are available direct water cleaning.
2. Keep water, liquids and humidity away from the machine, as moisture will damage the interior electronic components.
3. Store in cool, dry places. Avoid direct sunlight.

TECHNICAL DATA

VOLTAGE	DC 7.4 V
POWER	16W
WORKING TIME	15 Minutes
CHARGING	5V 1A

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

Malis B. Machoński Sp. K. declares that device described in this manual meets all the requirements of the:

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