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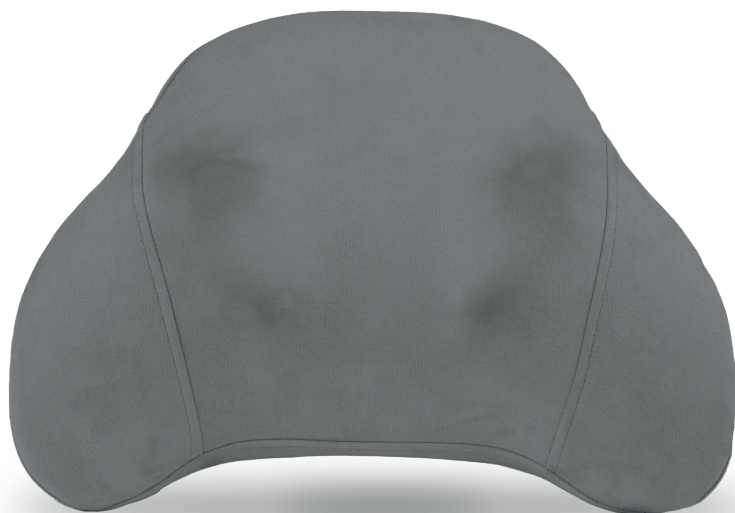
PL

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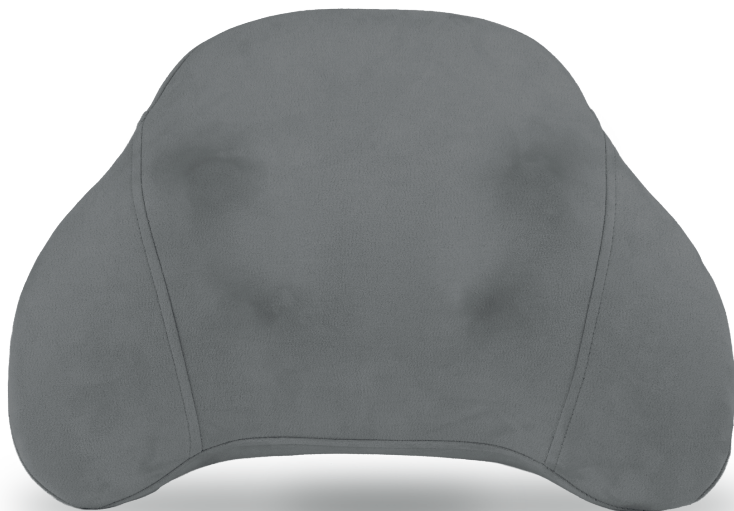
INSTRUCTION MANUAL EN

MASSAGE PILLOW
MAX-COMFY 1800 WIRELESS

Please read this user manual before using the device! Follow all safety instructions to avoid damage caused by improper use! Keep the user manual for future reference.

WARNING

To reduce the risk of fire, burns, electric shock, or other injuries, follow these instructions.

- Unplug the device from the power outlet after use.
- To avoid overheating, do not use the product under pillows, blankets, etc.
- Keep away from children.
- Keep away from water and liquids.
- Users with physical impairments or disabilities should use the product under supervision and after reading the instructions.
- Please use this product in accordance with the instructions provided.

DANGER

To reduce the risk of electric shock, follow these instructions:

- After use or before cleaning the product, unplug it from the power outlet.
- If the product falls into water, do not touch the wet product. You must unplug the power cord.
- Do not use the product while bathing in the shower or bathtub.
- Do not place the product in water or other liquids.

CONTRAINDICATIONS


1. Critical illness such as severe heart disease renal failure ascites due to liver cirrhosis severe systemic edema etc Prohibited.
2. Highly sensitive or severely allergic skin, contact skin disease patients are prohibited to use.
3. It is forbidden to use it where there are boils ulcers rashes and unexplained lumps on the body surface.
4. People with bleeding tendency such as late diabetes severe anemia leukemia, aplastic anemia and Thrombocytopenia etc. are strictly prohibited.
5. Low blood pressure, hypoglycemia, excessive weakness and nervous tension are strictly prohibited.
6. It is forbidden to use it if you are old, extremely weak emaciated or have lost elasticity of skin.
7. It is strictly prohibited to use in the painful parts or fracture parts of acute sprain and trauma.
8. Those who are too thirsty, too hungry, too drunk too tired or stay up too late are forbidden to use.
9. It is forbidden to use the lower abdomen during menstruation. It is

forbidden for pregnant women.

10. It is forbidden for patients with skin disorders or those who are not sensitive to heat.

11. Postoperative patients and those who install cardiac pacemakers are forbidden to use.

OPERATION METHOD

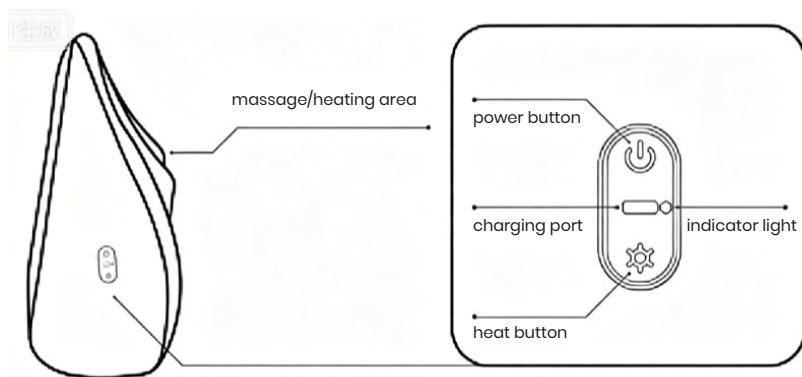Turning On and Off

Press and hold the Power button for 1 second to turn on the device.

The default settings after power-on are Massage Speed Level 2 and Heat Level 2. The green indicator light stays on.

To turn off the device, press and hold the Power button for 1 second.

Massage Speed Adjustment

The device has 2 massage speed levels. Briefly press the Power button to switch between the two massage speeds:

- **Level 1 – Low**
- **Level 2 – High**

Massage Direction

The massage direction changes automatically every 1 minute. Manual direction adjustment is not required.

Heating Function

The default heating setting after power-on is Heat Level 2. Briefly press the Heat button to switch between the two heating levels:

- **Level 1 – Low heat: indicator light is dim**
- **Level 2 – High heat: indicator light is fully illuminated**

Automatic Shut-Off

The device automatically turns off after 15 minutes of use.

CHARGING METHOD

1. This product cannot be used when charging
2. Plug the charging cable plug into the charging interface on the product.
3. When the product is in the working state, it is connected to the power supply for charging. The product starts the self-protection program and automatically stops working.
4. The indicator light is red when charging and green after charging.
5. Please unplug the power after charging.

Note: This product cannot be used during charging

Warning: the device has a heating surface, and people who are not sensitive to heat must pay attention when using it.

PRODUCT USE DIAGRAM

MAINTENANCE

Use a damp cloth to gently wipe the surface.

TECHNICAL DATA

POWER	25 W
VOLTAGE	5V
MATERIAL	Suede fabric
WIRELESS USAGE	7.4V / 1800mAh
RELAXING WARMTH	YES
SPEED CONTROL	YES
CHANGE OF DIRECTION OF THE HEADS	YES
LIGHT HEATING FUNCTION	YES
TIMER	15 minutes
ADJUSTABLE SEAT HANDLE	YES
INCLUDED	mains charger

WARRANTY

The seller is liable for defective goods for 2 years from the date of delivery.

All pictures and drawings are for reference only.

Disposal of electrical equipment

The symbol of the crossed-out wheeled waste container placed on the products informs that used electrical and electronic equipment (hereinafter also "waste equipment" or "electro-waste") may not be disposed of with municipal waste. By throwing e-waste into the container, you pose a threat to the environment. Used equipment may contain substances (e.g. lead, cadmium, chromium, bromine, mercury, freon) that are dangerous to human health and life and to the environment. By sorting and handing over used equipment for processing, recovery, recycling and disposal, you protect the environment against pollution and contamination, you contribute to reducing the use of natural resources and lowering the production costs of new equipment.

DECLARATION OF CONFORMITY

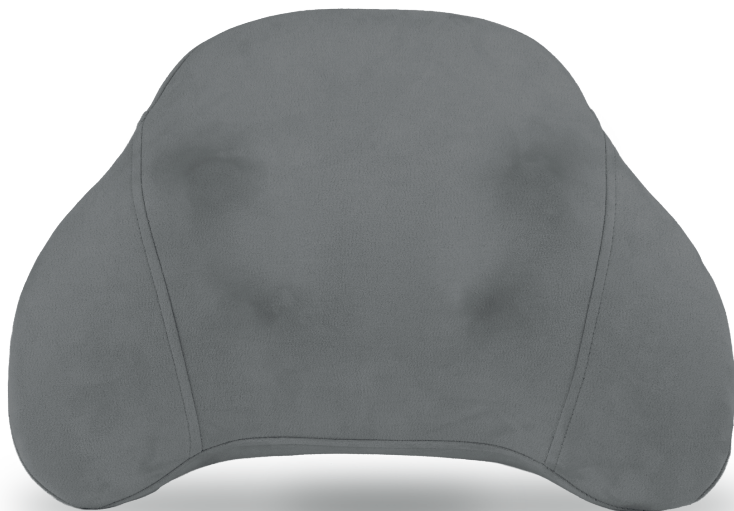
Malis B. Machoński Sp. K. declares that described in this manual Maltec meets all the requirements of the:

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